


Ultrassay Biotech Co., Ltd.

23th Floor, Building A, Future Tower, No. 711 Luzhou Avenue
 Hefei, Anhui, China
 E-mail: info@ultrassay.com Website: http://ultrassay.com



Disposable Cervical Sampler (A / B Type)

Disposable cervical sampler (A / B Type) is used for clinical collection of cervical cells for HPV examination and gynecological testing etc., It is designed according to the anatomical structure of women's cervix, the shape fits the cervical anatomy. the product is composed of PE brush head, PP rod and push tube, which is used for clinical sampling of gynecological cervical examination, with soft texture, no damage to the cervical mucosa, Using the push tube design, the brush head is pushed into the preservation medium bottle after sampling to avoid secondary contamination.
 Blister packing, EO sterile.
 10ml or 20ml preservation medium.

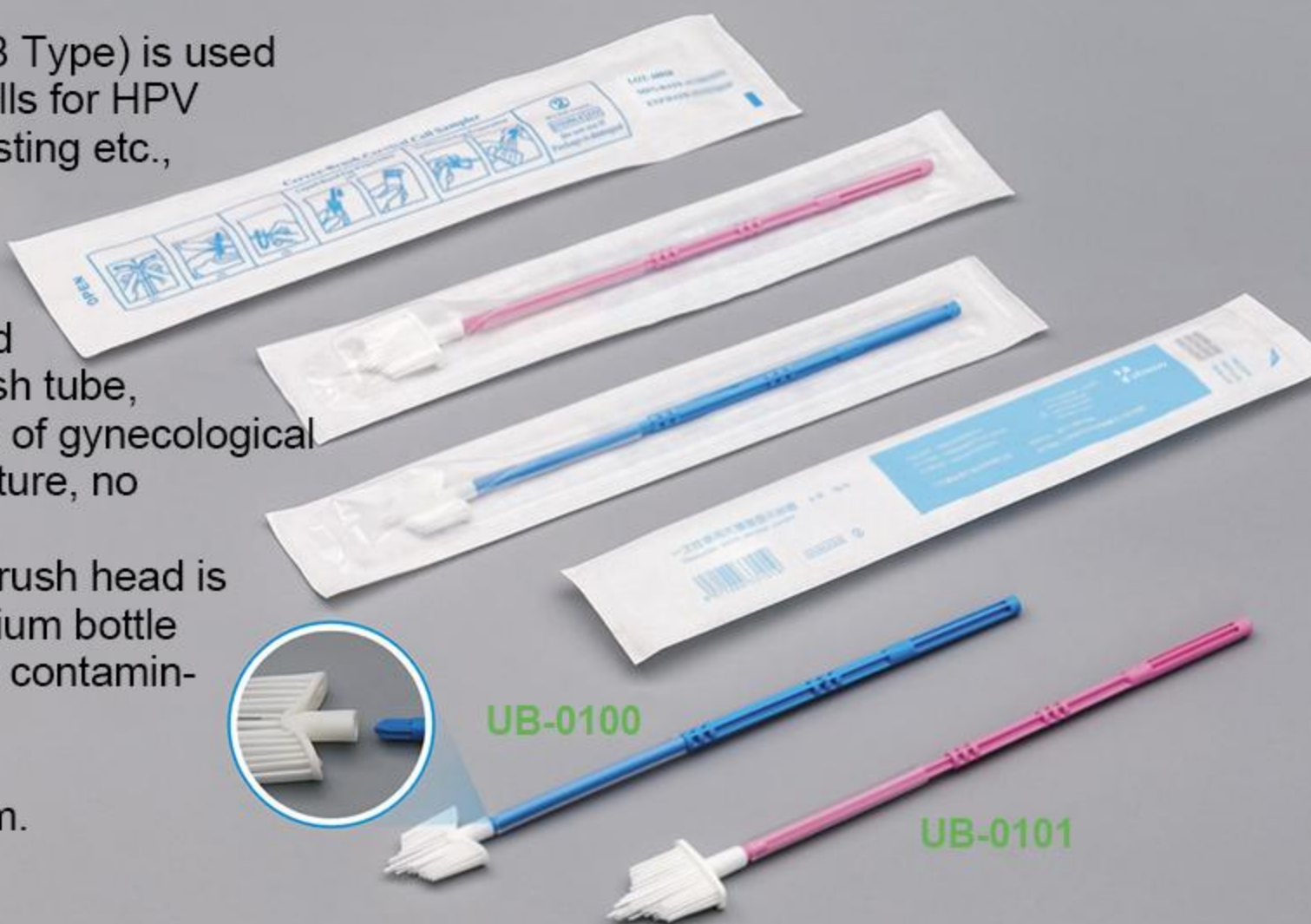

Ultrassay Biotech Co., Ltd.

Item No.	Product Name	Specification	Packing Size
UB-0100	Disposable Cervical Sampler	A Type	1000pcs/case
UB-0101	Disposable Cervical Sampler	B Type	1000pcs/case

Disposable Cervical Sampler (C Type)

Disposable cervical sampler type C is used for HPV examination, with soft texture, no damage to the cervical mucosa, It can collect more tissue samples, and simple operation and analysis. The breakpoint is designed at the rod of the type C sampler. and after sampling, the rod is broken to facilitate the brush head to be stored in the preservation medium for detection.
 Blister packing, EO sterile.



Ultrassay Biotech Co., Ltd.

Item No.	Product Name	Specification	Packing Size
UB-0102	Disposable Cervical Sampler	C Type (5.5cm breakpoint)	1000pcs/case
UB-0103	Disposable Cervical Sampler	C Type (7.0cm breakpoint)	1000pcs/case

Ultrassay Biotech Co., Ltd.

23th Floor, Building A, Future Tower, No. 711 Luzhou Avenue
 Hefei, Anhui, China
 E-mail: info@ultrassay.com Website: http://ultrassay.com



Self-sampling Kit

The self-sampling kit consists of a disposable cervical sampler type C and a preservation tube (containing the preservation medium).
 The HPV-2ML contains a 5ml cryovial with 2ml preservation medium;
 The HPV-3ML contains a 7ml cryovial with 3ml preservation medium.



Ultrassay Biotech Co., Ltd.

Item No.	Product Name	Specification	Packing Size
UB-0115	Self-sampling Kit	HPV-2ml	500pcs/case
UB-0116	Self-sampling Kit	HPV-3ml	500pcs/case

Universal Transport Medium (UTM)

Universal Transport Medium (Hank's) is used for virus culture without affecting virus activity and infection ability.
 Universal Transport Medium (PBS) is used to preserve the sample microenvironment and maintain the stability of the virus and other pathogens in the sample.
 Sample types: Nasopharyngeal/oropharyngeal swab and sputum.
 Store condition: Store in room temperature for 2 years.
 Sample storage: If the immediate sample handling is not required, keep the samples at -70°C for a long-time storage. Avoid freezing and thawing.



Ultrassay Biotech Co., Ltd.

Item No.	Product Name	Specification	Packing Size
UB-0121	UTM (Hank's)	5ml/7ml/10ml tube with 2ml/3ml medium	500pcs/case
UB-0122	UTM (PBS)	5ml/7ml/10ml tube with 2ml/3ml medium	500pcs/case



Disposable Swabs (Nylon Flocked)

Compared with traditional swabs, it has better sample adsorption and elution capabilities, suitable for secretions, viruses or cells samples collection. Made of medical grade ABS and nylon flocked.

Item No.	Product Name	Specification	Packing Size
UB-0111	Disposable Swabs		2000pcs/case

Disposable Swabs (Sponge)

It is made of medical grade ABS and sponge with 5cm breakpoint, which is used for the collection of various secretions from the human body.

Item No.	Product Name	Specification	Packing Size
UB-0112	Disposable Swabs	5cm breakpoint	2000pcs/case

Disposable Cervical Sampler (Self-sampling Kit)



UB-0118



UB-0119



UB-0120

The disposable Cervical Cell Sampler (Self-sampling Kit) is used for women to collect cervical cell samples at home. The disposable Cervical Cell Sampler (Self-sampling Kit) is composed of protected shell, rotatory rod, brush head and cap head. The samplers are divided into A, B, and C types according to the style of brush heads. Type A is made of medical grade nylon flocked, Type B is made of medical grade PE, Type C is made of medical grade sponge. EO Sterile.

Item No.	Product Name	Specification	Packing Size
UB-0118	Disposable Cervical Sampler	A Type	800pcs/case
UB-0119	Disposable Cervical Sampler	B Type	800pcs/case
UB-0120	Disposable Cervical Sampler	C Type	800pcs/case