


Available rotors



Swing bucket smear rotor  
Max. Capacity: 6/8pcs  
Max. Speed: 4000rpm  
Max. RCF: 2170xg  
Order No.: CYTOR-1/2



Bucket  
Order No.: CYTOF-1

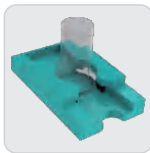


Plastic holder  
Order No.: CYTOR-2



Cell smear rotor  
Max. Capacity: 6pcs  
Max. Speed: 3000rpm  
Max. RCF: 730xg  
Order No.: CYTRO-4



Smear Clip  
Order No.: CYTOF-3

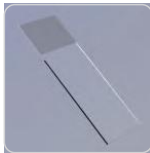

Glass slide  
Order No.: CYTOF-4



Cell smear rotor  
Max. Capacity: 12pcs  
Max. Speed: 3000rpm  
Max. RCF: 890xg  
Order No.: CYTRO-5



Filter paper  
Order No.: CYTOF-5



Swing bucket smear rotor  
Max. Capacity: 16pcs  
Max. Speed: 2000rpm  
Max. RCF: 940xg  
Order No.: CYTRO-3

Technical specification

Model	Cytoprep-1	Cytoprep-2
Max. Speed	4000rpm	2000rpm
Max. RCF	2170×g	842×g
Available Rotor	6 or 8 pcs	12 or 16 pcs
Time Range	0~99mins69s	
RPM/RCF Convert	Yes	
Noise	≤ 55 dB	
Acc/Dec	10 Kinds	
Speed Accuracy	±20r/min	
Display	LCD	Touch Screen
Memory	10 groups	100 groups
Power Supply	AC 110~220V, 50/60Hz	
Overall Dimension	528×370×280mm	690×500×390mm
Net Weight	21kg	41kg
Testing Standard	EN 61010-1, EN 61010-2-020, EN 61326-1, EN 61010-3-2/A2	



**Ultrassay BioTech Co., Ltd.**  
RM3018 Poly D9, No. 999 Luzhou Avenue, Baohe District, Hefei, Anhui, China  
Tel: +86-551-62881663  
E-mail: info@ultrassay.com  
Website: https://ultrassay.com

Cytoprep  
Cell smear Centrifuge

Applications

ULTRASSAY Cytoprep cell smear centrifuge can be used for the liquid handing in fields of cytology, hematology and microbiology. Such as cerebrospinal fluid, pleural fluid, peritoneal fluid, synovial fluid, urine, oral cavity-cleaning, stomach cleaning, exudate, percolate, phlegm. It is mainly used for checks with small cell capacity. Our clip rotor for cytology is compatible with popular cytocentrifuge accessories currently on the market.



Cytoprep-1

Features

- **Convenient Operation**
  - (1) Special slides preparation device, one-step rapid handling, easy slides unloading, time saving and labor saving;
  - (2) One-time automatic slides preparation of 1~16 pcs, which takes only 3 minutes;
  - (3) 6, 8, 12 or 16 through-put available for selection.
- **High Compatibility**
  - (1) The cell preservation solution is suitable for a variety of cytological specimens;
  - (2) Rapid staining program;
  - (3) Papanicolaou and hematoxylin-eosin staining available.
- **High Quality Slides Smear**
  - (1) 1.8cm diameter slides smear, convenient to read;
  - (2) Multiple slides of the same quality can be made repeatedly from the same specimen;
  - (3) The background of the slides smear is clear, and the cells are abundant and evenly distributed.
- **Main Consumables**

Liquid-based thin-layer cells slide preparation instrument, cells preservation solution, cells filter, staining solution, etc.



Cytoprep-2